




<b>Product</b> <b>Item number</b>	CareSystem support arm, fixed height R3018000
<b>Colour</b>  <b>Product description</b>	<div> <div>White</div>  </div> <p>Support arm without spring counter balance, 850 mm, fixed height, folds against the wall. Incl. mounting kit for brick wall. Discontinued as per 31 Dec 2015</p> <p>Support arm without spring counter balance, 850 mm, fixed height, folds against the wall. Incl. mounting kit for brick wall. Discontinued as per 31 Dec 2015</p>
<b>Adjustability</b> <b>Product capacity</b>	Fixed mounting Max. vertical load: 135 kg.
<b>Materials</b>  <b>Ergonomy</b>  <b>Features</b>	<p>Basic material: 20 micrometre anodized aluminium</p> <p>Surface material: Two-component plastic (PVC)</p> <p>End cap/cover cap: Arnite</p> <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p> <p>Folds against the wall.</p> <p>The underside is shaped for the mounting of various accessories.</p>

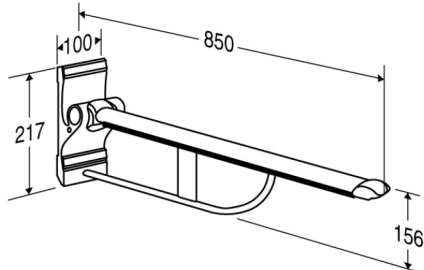

# PRESSALIT®

<b>Product</b>	CareSystem support arm, fixed height
<b>Item number</b>	R3018000
<b>Cleaning</b>	Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

PRESSALIT®