




Product Item number	CareSystem support arm with spring fixed height RF016000
Colour Product description	<div> <div>White</div>  </div> <p>Support arm with spring counter balance, 600 mm, fixed height folds against the wall. Incl. mounting kit for brick wall. Discontinued as per 31 Dec 2015</p> <p>Support arm with spring counter balance, 600 mm, fixed height folds against the wall. Incl. mounting kit for brick wall. Discontinued as per 31 Dec 2015</p>
Adjustability Product capacity	Fixed mounting Max. vertical load: 135 kg.
Materials Ergonomy Features	<p>Basic material: 20 micrometre anodized aluminium</p> <p>Surface material: Two-component plastic (PVC)</p> <p>End cap/cover cap: Arnite</p> <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>The surface material is soft and slip-resistant, and has a warm feel.</p> <p>The spring ensures that the support does not "fall" down from its vertical position, and it makes raising easier. A counterbalanced support arm can be raised and lowered with a sliding movement or in easy stages. Folds against the wall.</p> <p>The underside is shaped for the mounting of various accessories.</p>

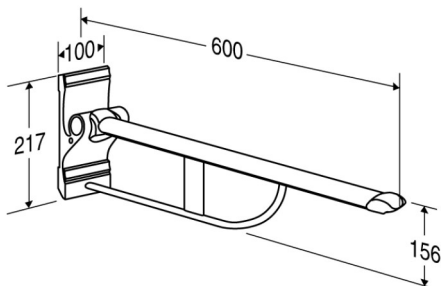

PRESSALIT®

Product	CareSystem support arm with spring fixed height
Item number	RF016000
Cleaning	Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®