




Product Item number	PLUS backrest, fixed height R6404112
Colour	Anthracite grey
Product description	<p>VVS: 777530212 NRF: 6317034 EAN: 5708590365398 HMI: 124215</p>  <p>PLUS backrests are offered in a number of colours and versions, they are available as both fixed height and adjustable. They are suitable for wall-mounting behind a toilet without cistern or for mounting on one of Pressalit's toilet lifters.</p> <p>PLUS backrest, 385 x 209 mm, soft feel. Fixed height, for mounting directly on the wall behind a wall-mounted elongated toilet without cistern. Materials: Aluminium and polyurethane foam with a plywood core. Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p>
Adjustability Product capacity	<p>Fixed mounting</p> <p>Maximum load: 50 kg. Maximum user weight: 100 kg. When the backrest is loaded with 50 kg, each wall fastener will be carrying 118 kg. <i>Backrests need to be able to support the weight of the person who is leaning against it.</i> <i>The backrest is not intended to carry the full body weight of a person, and the test methods applicable to backrests reflect this: ISO 17966 prescribes a static test at 50% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 50 kg, it is – according to ISO 17966 – the equivalent of a user weighing 100 kg leaning on the backrest.</i></p>
Materials	<p>The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer.</p>
Ergonomy	<p>Soft and slip-resistant material that has a warm feel.</p>
Features	<p>For toilet without cistern</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>

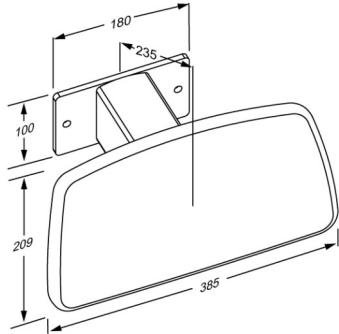

PRESSALIT®

Product	PLUS backrest, fixed height
Item number	R6404112
Cleaning	<p>PLUS backrest.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
Declaration of conformity	<p>CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices.</p> <p>Expected service life: 10 years.</p> <p>Contraindications for use: It is not recommended that people who need to lie down while being changed use the backrest for support.</p>

pressalit.com

PRESSALIT®